

■ COMPACT SOLUTION FOR BIOHAZARD WASTE DISPOSAL

Integrated Sterilizer & Shredder (ISS)

The ISS performs both shredding and waste steam sterilization in a single vessel. The vessel is fitted with a motor-driven shaft, with powerful shredding/crushing blades which reduce the size and volume of the waste.



- Easy to Operate – No need for special technician qualification.
- Environmentally Sound - Shredded waste is reduced to as little as 1/5 its original volume, without emitting harmful substances.
- Cost-effective - Inexpensive operation and maintenance.
- Totally Safe - Automatic locking door prohibits unauthorized interruption.
- Efficient - A single unit can serve any middle size hospital, clinic or laboratory.
- Compact - Room of only 3 x 4 meters is necessary.
- Easily and quickly installed – Within one day.



ISS AC-575



PharmaTherm Limited
Ireland & UK Distributor of Celitron Autoclaves
Address: Carrignavar, Cork
Telephone: 00353 (0)21 4884837
Web: www.pharmatherm.ie
E-mail: sales@pharmatherm.ie



Celitron Medical Technologies Kft.
Address: H-2600 Vác, Szent László str. 36. Hungary
Telephone: +36 27 412-610, Fax: +36 27 412-611
Web: www.celitron.com
E-mail: info@celitron.com



high-tech in sterilization
Professional Sterilizers and Medical Waste Solutions




Celitron's compact medical hazardous waste solution, the Integrated Sterilizer & Shredder, is a large steam sterilizer with an integrated shredder, designed for on-site conversion of biohazard waste in hospitals and clinics, complying the EU and WHO recommendations.

The stainless steel vessel is fitted with a motor-driven shaft, with powerful shredding / crushing blades that can rotate in two directions to reduce the size and volume of the waste down to 20%. The 5.5kW motor is sufficient to rotate the shaft with an RPM of 400-1700 for various operations. The blades are mounted on the shaft and are designed to shred waste such as sharps, dialyzers, syringes, papers, cloth, plastic and glass. Shredding is important as it enhances the steam penetration therefore improves the overall sterilization results.

The steam is internally produced by a steam generator, supplied by water purification and draining system as standard accessory. The vessel is constructed with internal sprinklers for automatic cleaning.

The entire process is automatic, including the opening and closing the door, turning and rotating of the vessel and the sequences of shredding and sterilization. The total cycle time can take as fast as 25 minutes and between 15 to 25 kg of waste can be processed.

After treated in the ISS, the waste is sterile. The liquid components of the waste are steamed out of the vessel, re-condensed and drained to a municipal sewer. As the waste is dehydrated, there is no risk of contaminated waste water. The waste is rendered, fragmented, non-toxic, largely solid and dry and therefore safe to be disposed as regular municipal waste.

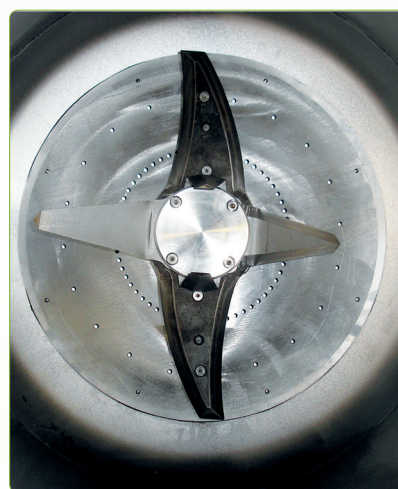

SPECIFICATIONS:

Model ISS	ISS AC-575
Chamber dimensions Inner dia. x Depth	500 x 800 mm
Chamber volume	150 l
Approximate weight	880 kg
External dimensions W x H x D	1290 x 2150 x 2039 mm
Chamber door	Automatic door with advanced safety features
Available sterilization cycles	Waste 134 °C / Textiles no Cut 134 °C / Special Waste 134 °C
Available test and other cycles	Dynamic test / Cleaning cycle
Steam source	Built-in steam generator (with a possibility of external steam connection)
Steam pressure (relative)	3 bar (43,5 psi)
Compressed air	6.0 – 8.0 bar (87 -116 psi)
Water source	Filtered Tap Water
Water pressure	1.0-6.0 bar (14.5-87 psi)
Power supply*	3-Ph. 380-400 V, 50/60 Hz
Power with steam generator	25 kW

* Adjustable to different voltage systems